



Central Display Unit

Nurse - Call Systems



Emergency Toilet Unit



Door Display Unit


Medisystems Electronic Nurse-Call Systems employ the latest technologies to deliver a number of useful features : Call Acknowledgment, Call Reminder, Nurse Presence Registration, Emergency Alerts, Nurse Help Requests, Instrument Alarm Relay, I-V Drip Alert, Hospital-Wide Code-Blue Alert, SMS Alert, Patient-Nurse Intercom, Corridor Displays, Call Transfer, Attendant Call, Configured Handsets offering linkage with Lighting, Audio and TV facilities, System Diagnostics and Nurse-Call Response Monitoring. These systems are available in several models and can be configured for every kind of hospital.



Bed Unit



Patient Handset

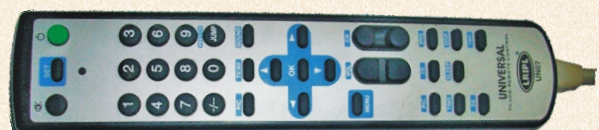

Patient Universal Handset

Bed - Head Panels



ICU Bed - Head Panel

Medisystems Bed-Head Panels converge every needed hospital utility at the patient's bed. Fully pre-fabricated aluminium panels, ensure organised convergence of Medical Gas Outlets, Electrical Mains, Standby / UPS power points, Nurse-Call Units and Bio-Signal, Data & Intercom outlets. International Standard Medical Rails in Stainless Steel, along with a full range of rail mountable accessories, enable every facility to be neatly made available at the patient's bed. Panels can be custom configured for every type of ICU, Patient Room, General Ward or Examination Room.



Vertical Panels



Horizontal Running Panels



NICU & Paediatric Panels



Concealed Panels

Routine add-ons include patient-bed lamps available in two and three feet lengths along with panel mounted switching.



Patient - Bed Lamp

Aesthetically designed Bed-Head Panel and Patient-Bed Lamp combinations for Deluxe Rooms, Private Rooms and Suites. Suitable versions are also available for Semi-Private Rooms and General Wards.



Rotating Pendant



Rigid Pendant

OT Pendants

Medisystems OT Pendants are fabricated in pure Stainless Steel and Aluminium. They are available in Rigid or Rotating versions in Shell or Shelf configurations. All electricals are protected and wired in sealed and fail-safe mode. Medical gas outlet arrangements are customised to user's choice.



Single Arm Pendant

OPD Patient - Call Systems



Main Display Unit



Desk Unit

Medisystems OPD Patient-Call Systems address the daily chaos in OPD waiting areas in hospitals, everywhere. Installation of these systems brings a huge improvement in the management of OPDs. The system can address upto 32 examination rooms, consulting rooms or counter windows. The main display can also be slaved for simultaneous visibility in different waiting areas. No

central operator is needed. The physician directly keys in the patient / case number and the door display indicator helps guide the patient into the room.



Door Display

CR MEDISYSTEMS PVT. LTD.

335 Pragati Estate, N.M. Joshi Marg, Lower Parel, Mumbai – 400011, INDIA. Tel: 2300 4930, 23094416
Fax: 91-22-2306 1903 • E-mail : medisystems2@gmail.com • Website : www.medisystems.in