

MEDISYSTEMS / 3200 – MEDICAL RAIL & ACCESSORIES



Medisystems offers a range of medical rail mountable accessories. Each of these accessories is fixed on one or more sliders which can be smoothly moved on the rail and fixed at any desired location along the length of the rail.

1. **Case Sheet Holders** > Medical case sheet holders come in two sizes, small and large. The large one is sized to hold X-Ray and Sonography films as well. Case sheet holders are available in the standard version – powder coated aluminium, or in 304 grade stainless steel.
2. **Bowl Holders** > These are ring holders assembled to hold a standard 6” bowl or a larger 8” bowl. The ring holders and bowls are made of 304 grade stainless steel.
3. **Utility Basket Holders** > These are handy stainless steel wire baskets for general purpose miscellaneous items. They come in two sizes small and large.
4. **Blood pressure Instrument holders** > Mercury or Dial Gauge Instruments can be mounted on the medical rail slider for convenient location and immediate availability behind the patient’s bed.
5. **I-V Pole Assembly** > A height adjustable stainless steel pole is mounted on a slider assembly. The pole is crowned with a hook arrangement for up to two I-V bags. The availability of an I-V pole on a medical rail helps reduce floor space clutter around the patient’s bed.
6. **I-V Pole Assembly w/Infusion Pump Stand** > This assembly includes an Infusion pump Stand mounted on the vertical limb of the I-V pole. The pole diameter is sized to mount all international brands of syringe and infusion pumps. A stack of four to six pumps can be accommodated.
7. **Suction Bottle Holder** > This is a simple clamp to enable mounting of the suction bottle below the vacuum terminal when medical gas outlets have been brought out through the bed head panel or the pendant.
8. **Examination Lamp Assembly** > An articulated examination lamp is mounted on a slider for convenient location behind the patients bed. The examination lamp is fitted with a self- reflective led lamp. When not in use, it can be swung out of the way but is always available for use when needed.
9. **Rail Slider Clamp** > This is a bare slider which fits and slides on the medical rail. It can be used to mount any object or accessory of convenience onto the rail.

10. **Syringe Pump Holder** > Designed for stacking 3-4 syringe pumps vertically on a single rail clamp.
11. **Patient Monitor Stands** > Designed to securely hold patient monitors in fixed, vertical height adjustment and wall angular positions. Come with specific clamps for securely holding any international brand of patient monitor on the panel or wall. Can also come with a tray for placing monitors without fixing. Monitor cable basket is included.

* * * * *

MEDISYSTEMS / 3200 – MEDICAL RAIL ACCESSORIES

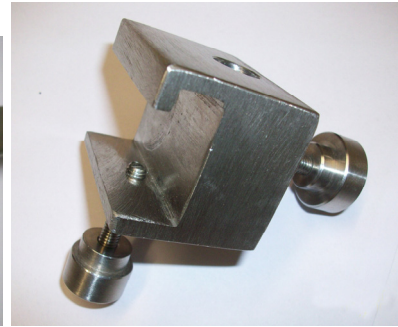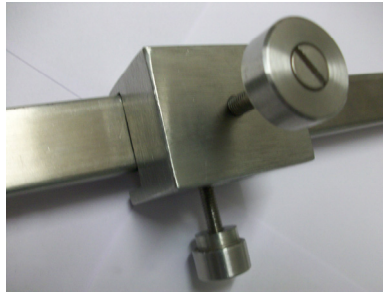


Medical Rail S.S. 304(25mm*10mm) Rail Slider Clamp S.S.304 with(13mm φ)



Case Sheet Holder std. 6" Bowl Holder with clamp

I.V. Pole Assembly w/clamp

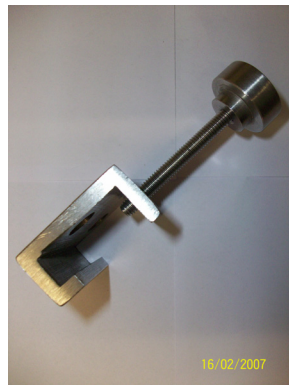

IV slider clamp (13mm φ)

Rail Slider Clamp Utility Basket Holder

Suction Bottle Holder



LED Examination Lamp



***Syringe Pump
Holder (25mm ϕ)***



***Syringe Pump/
Monitor Stand Holder (25mm ϕ)***



Bed-Head Panel with Rail & Accessories